




REHASENSE

SOFTGUARD

PRESSURE REDUCING MATTRESSES

HIBERMATE

The Hibernite pressure reducing mattress powered by Mammoth Technology has been designed not only to offer optimum pressure reducing qualities but has also been clinically proven to enhance sleep efficiency.

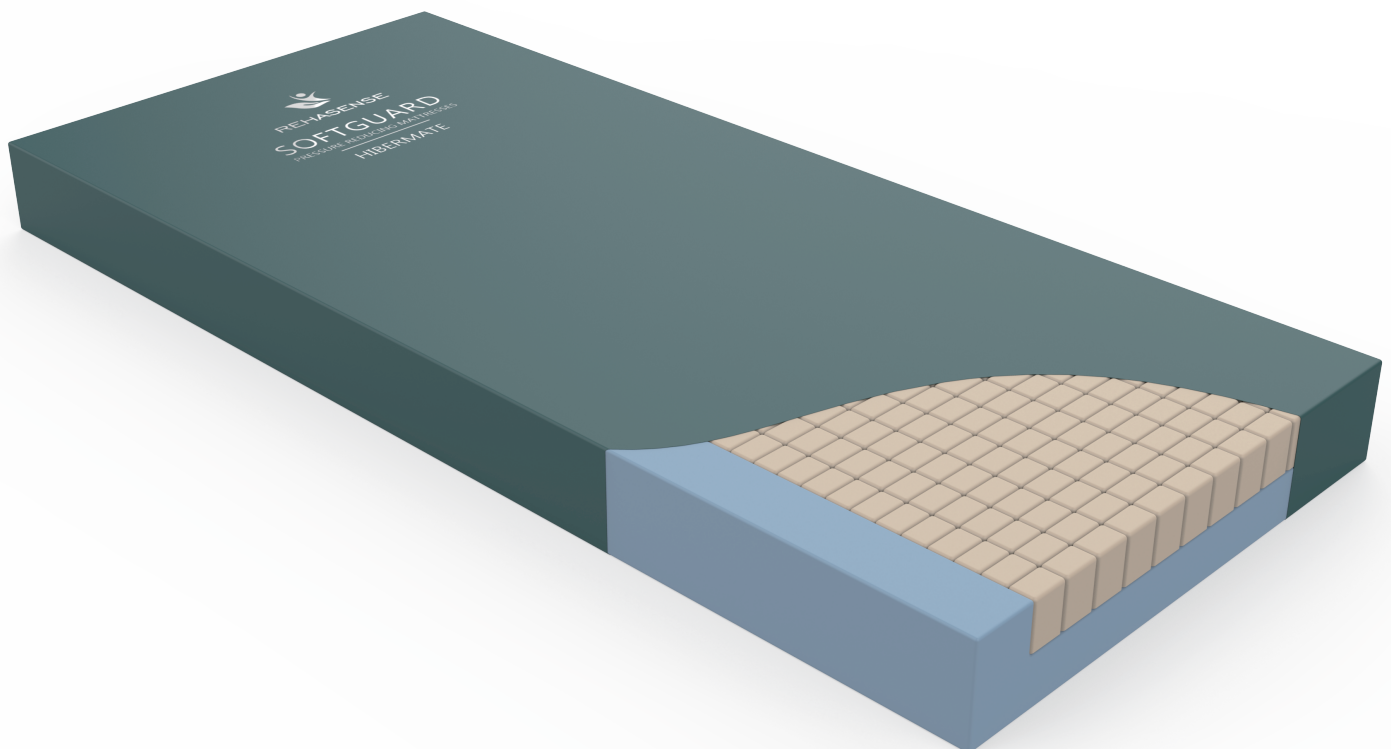

SOFTGUARD

PRESSURE REDUCING MATTRESSES

HIBERMATE



REHASENSE

KEY FEATURES

The Hibermate pressure reducing mattress features the Mammoth® V-Max® castellated cells which combine a specific foam density with a nominal hardness that optimises pressure reduction over high risk areas. The subsequent improved circulation through the soft tissues not only reduces the risk and management of pressure ulcers but has demonstrated to improve the quality of sleep.



- The high specification foam contributes to an effective pressure ulcer prevention care plan
- Internal air channels in the foam's surface help to reduce the build-up of heat and moisture on the skin
- Smooth multi-stretch cover allows the patient to glide over the surface of the mattress
- The waterproof, vapour permeable cover can be easily wiped clean
- High frequency welded seams reduce the risk of fluid ingress
- Full length welded flap cover can be removed and laundered at 95°C

KEY FACTORS

Tested and developed under controlled scientific conditions, 3 key sleep factors were demonstrated to improve:

Sleep efficiency

Participants improved sleep efficiency by a full clinical banding, moving up from 'poor sleep', to 'good sleep' thanks to a 7% improvement.

Sleep onset latency

Participants fell asleep quicker showing a 28% improvement.

Sleep comfort & enjoyment

Improved by 21%. Improved sleep efficiency means better healing times, shorter hospital stays, and a reduction in wound management costs.

TECHNICAL INFORMATION

Product Code:

Dimensions:

Single 198 x 88 x 15.6 cm

King Single: 198 x 105 x 15.6 cm

Maximum Weight Limit:

36 stone / 229 kg

Risk Factor:

High Risk

Warranty:

3 years cover / 2 years foam

Mattress Weight:


14 kg

Safety Standards:

BS EN 597-1 Cig

BS EN 597-2 Match

BS7177: 2008+A1 Crib Source 5

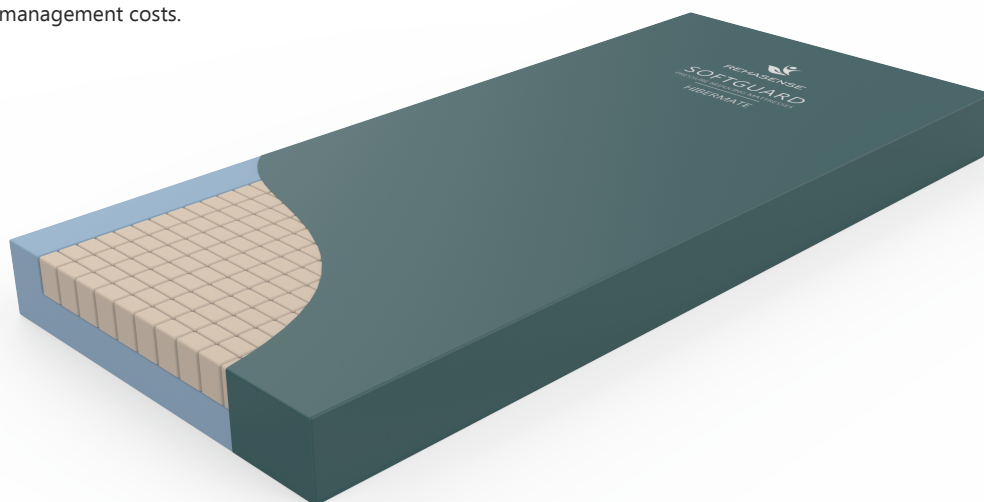

The **Hibermate** is CE-marked and fulfils the requirements of EU directive 93/42 EEC for medical devices.

Copyright © Rehasense® 2019